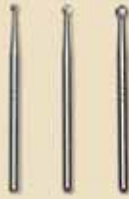


Rotary

Evian Perio Burs

Surgical Length FG 25mm for Periodontal Debridement

- Blunt, tissue-kind heads
- Removes granulation tissue efficiently
- Removes deposits from root surfaces
- Easier technique than the tedious task of hand scaling



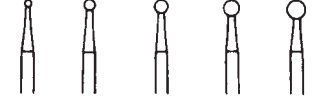
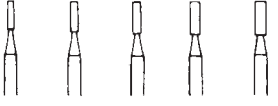

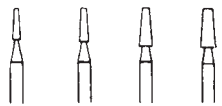
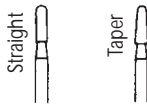
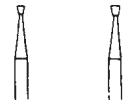


Package of six

9052041	2
9052042	4
9052043	6
9052040	Assorted (2 of each size)

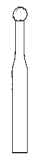





Surgical Length Carbide Burs (FG Surgical 25mm) *Made in Germany*

1 per package

 Round <i>Up to #2 for single surface cavities, #4 for anterior interproximals. In endodontics, ideal for opening pulpal chamber.</i>		 Cross Cut Fissure Straight <i>To produce parallel walls and flat floors in preparations</i>					
 1 2 4 6 8 008 010 014 018 023 SFG SFG SFG SFG SFG		 556 557 558 559 560 008 010 012 014 016 SFG SFG SFG SFG SFG					
 Cross Cut Fissure Taper <i>To avoid undercuts and create slightly divergent cavity walls</i>		Breakthru Burs <i>To produce undercuts</i>	Inverted Cone <i>For large occlusal restorations to insure rounded internal line angles</i>				
 700 701 702 703 010 012 016 021 SFG SFG SFG SFG		 Straight Taper 1558 402 014 016 SFG SFG	 34 35 008 010 SFG SFG				
9050601	1	9050607	34	9050618	556	9050625	700
9050602	2	9050608	35	9050619	557	9050626	701
9050603	4	9050612	402	9050620	558	9050627	702
9050604	6	9050615	1558	9050621	559	9050628	703
9050605	8			9050622	560		

Regular Steel Burs (HP 44.5mm, RA 22mm) *Made in Germany*

10 per package

							
	Round	Inv. Cone	CC Sq. Fis.	CC Tap. Fis.	LH Round	Wheel	
HP	1/2, 2, 4, 6, 8	—	558, 562	700	10, 12	—	
9052051	1/2	9052057	6	9052083	562	9052063	12
9052053	2	9052059	8	9052086	700		
9052055	4	9052079	558	9052061	10		
RA	1/2, 2, 3, 4, 5, 6, 7, 8	35, 37	—	—	—	14	
9052001	1/2	9052005	4	9052008	7	9052022	37
9052003	2	9052006	5	9052009	8	9052016	14
9052004	3	9052007	6	9052020	35		