

# Gynecology

["LOOP" Surgery Products](#)

[Medicaments](#)

[Pessary Products](#)

[Biopsy Forceps](#)


[Specula](#)

[Uterine Dilators](#)

[Cervical Dilators](#)

[Uterine Sounds](#)

[Instruments](#)

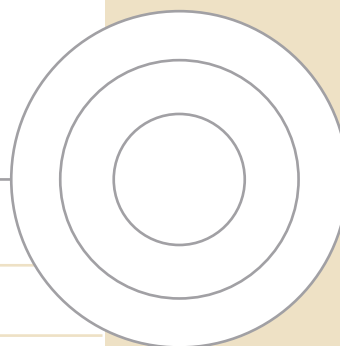



*Designed with  
the alternate site  
in mind!*

### Model 2001e Electrosurgical Unit

Full-Featured Electrocautery  
for the Office

- 120 watts Maximum Power Output  
*Perform office-based procedures with confidence.*
- Five Distinct Output Modes  
*Cut, Coagulation, Blend, Fulguration and Bipolar outputs provide unsurpassed versatility.*
- Isolated Output Transformers and Patient Circuit Monitor  
*Ensure safety to both patient and staff.*
- Flush Membrane Front Panel  
*No knobs or switches to adjust, easy to clean and disinfect.*
- Digital Diagnostic Error Codes  
*Readily identify and protect against potential fault conditions.*
- Hand Control with Optional Foot Control  
*Fully adaptable to almost any procedure or electrode.*
- Meets the Following Safety and Performance Standards:  
*CSA C22.2 NO 601.1-M90, UL 2601-1-UL, IEC 60601-2-2 (1998-90)*

9006501 Model 2001e Electrocautery Unit (110/240VAC)  
(Includes: hospital-grade power cords (110V and 220V) disposable hand-control pencil, disposable patient pad, clamp cord assembly and manual)

**SPECIFICATIONS:**

Maximum Power Outputs by Operating Mode:  
 Pure Cut = 120 Watts @ 500Ω, 357kHz ± 50 kHz  
 Blend = 90 Watts @ 800Ω, 357kHz ± 50 kHz  
 Coagulation = 80 watts @ 1000Ω, 475kHz ± 19 kHz  
 Fulguration = 40 watts @ 1000Ω, 410kHz ± 50 kHz  
 Bipolar = 30 watts @ 200Ω, 535kHz ± 29 kHz

AC Input Current: 110/240VAC  
 Mains line frequency range (nominal): 50-60Hz  
 Power consumption: 540VA

Dimensions:  
 Size – 10 1/4”(W) x 12”(L) x 6”(H)  
 (26cm x 30.5cm x 15.2cm)  
 Weight – 14 lbs (6.5kg)

Two (2) year warranty



## Evacuator/Air Filtration System

The Premier Evacuator/Air Filtration System with its distinctly quiet operation (60 decibels at 1.0 meter) is suitable for electro-surgical or laser procedures. ‘Responsible’ filter pack conveniently and safely collects plume (smoke and particulate matter to 0.2 microns at 99.999% efficacy<sup>1</sup>) and neutralizes odors. Each filter pack comes with 3/8” tubing to attach to specula fitted with a smoke-removal tube.

1. Filter rating is for 15 procedures or 45 minutes of accumulated filtration time, whichever comes first.

- 9006300 Premier Evacuator (110V)  
Includes One (1) Filter

---

- 9006298 Premier Evacuator (220V)

---

- 9006290 Replacement Filter Packs  
(5 ea/case)

---

- 9006292 Reusable Smoke Tube  
(for accessible sites)

---

- 9006296 Evacuator Footswitch  
(not shown)

---

### SPECIFICATIONS:

AC Input Current:

- 9006300 - 110/120VAC, 60Hz, 15 amps
- 9006298 - 220/240VAC, 50Hz, 8 amps

Dimensions:

- Size - 10”(W) x 10”(L) x 14”(H)  
(25.4cm x 25.4cm x 35.6cm)
- Weight - 20lbs. (9.1kg)

Two (2) year warranty

## “LOOP” Surgery Products

### ESU/Evacuator Procedure Cart

- Constructed of durable, non-conductive ABS plastic
- 300 lb. (135kg) capacity
- Three large, spacious shelves
- Transport handle and quiet, directional casters included
- Snap-together assembly (no tools required)

9006560 Premier ESU/Evacuator Procedure Cart

Dimensions:

Shelf Size - 16" x 24"  
 (40.64cm x 60.96cm)  
 Height - 34 7/8" (88.6cm)  
 Weight - 32 lbs. (14.4kg)  
 Capacity - 300 lbs. (135kg)

*Optional procedure cart pulls the whole package together!*



## Hand-Control Pencil



Activates generator's cut and coagulation output modes. Fitted with a 10' cord, accepts standard shaft size (3/32") electrodes. (Compatible with other manufacturers' generators.)

9006505 Hand-Control Pencil, Disposable, 5/pack

9006515 Hand-Control Pencil, Reusable, one each

## Foot-Control Pencil



For use with Foot Controller (9006525). Ten foot cord comes with standard banana plug (Bovie plug adapter needed for use with 2001e). Accepts standard shaft size (3/32") electrodes. (Also can be used with other manufacturers' generators.)

9006530 Foot-Control Pencil, Disposable, 5/pack

9006585 Bovie Plug Adapter, one each

## Dual-Foot Controller



Dual operating mode conveniently selects either cut or coagulation operating modes. Standard color-coded pedals and fluid resistant cover (not shown) included. Also for use with bipolar accessories (see page 28).

9006525 Dual-Foot Controller

*(For use with the Premier 2001e generator only.)*

## "LOOP" Surgery Products

### Disposable Dispersive Patient Pads



9006540 Dispersive Pad<sup>2</sup>, Single Cell, Disposable, 5/pack



9006550 Dispersive Pad<sup>3</sup>, Dual Cell, Disposable, 5/pack  
For use with pad monitoring systems.

9006535 Replacement Clamp Cord Assembly (not shown)  
For connecting unrecorded patient plates to the Premier 2001e generator.

### PolyPlate™

Reusable Dispersive Patient Plate



**Innovative** - Utilizes advanced material technology.

**Safe** - Eliminates the risk of skin rashes, allergic reactions and possible lesions that can occur with the use of adhesive pads or pads requiring gels.

**Economical** - Saves costs, easy to clean and reduces medical waste.

General usage guidelines are printed on the non-patient contact side of the plate. Package insert provides detailed product instructions.

For use with all electrosurgery generators without a dispersive plate monitoring system. Compatible with Premier and non-Premier generators.

9006545 Reusable PolyPlate™, Single Cell

2. All Premier dispersive pads are unrecorded.  
3. Not for use with 2001e.

## “LOOP” Surgery Products





### TruCone™ Rotational Cone Biopsy Electrode

- Easy to Use  
*No awkward hand movements required.*
- Physician Assured  
*Stabilizing prongs increase procedure confidence.*
- Improved Diagnosis  
*Produces a clean, smooth cut with minimum thermal artifact. Sample conveniently secured.*
- Patient Friendly, Safe  
*Insulated components ensure safety to the patient, physician and staff.*
- Eight (8) Different Sizes  
*The largest number of sizes available.*

Packaged five per box. Each TruCone is packaged with instructional pathological guidelines.

	Diameter	Depth
9006110	18mm	12mm
9006115	18mm	15mm
9006125	24mm	12mm
9006130	24mm	15mm
9006140	30mm	12mm
9006145	30mm	15mm
9006165	30mm	9mm
9006170*	12mm	21mm

9006175 Assorted Pack (3)  
1 each: 18mm x 12mm, 18mm x 15mm, 24mm x 15mm

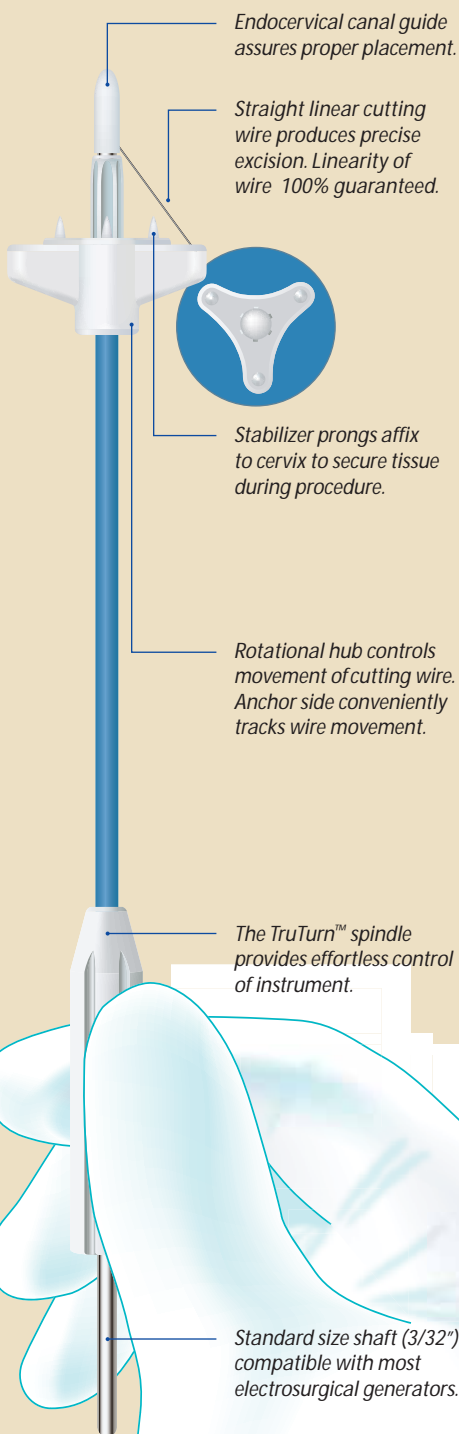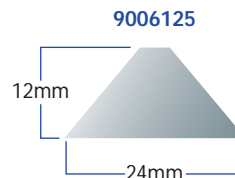

## Obtain Precise Cervical Cone Biopsies – Effortlessly! *It is as easy as...*

1. With the electro-surgical generator in the active mode, insert the appropriately sized TruCone™ into the endocervical canal. (Use of a foot-controlled generator recommended.)
2. After the stabilizing prongs have secured the cervix, gently rotate the shaft spindle for one complete rotation. Hub serves as a rotation guide.
3. Deactivate generator and remove TruCone instrument containing secured sample from vaginal cavity.

\* Stabilizer prongs not available on this size. / US Patent No. 5,676,663

## "LOOP" Electrodes

All Premier Electrodes are made with tungsten wire and fitted with a standard  $\frac{3}{32}$ " diameter, 5" long shaft. Compatible with most electrosurgical generators on the market.


### Reusable "LOOP" Electrodes

5 to 6 uses with proper care. Available individually.



### Disposable "LOOP" Electrodes

Sterile, single-use electrodes. Available in 5-packs.

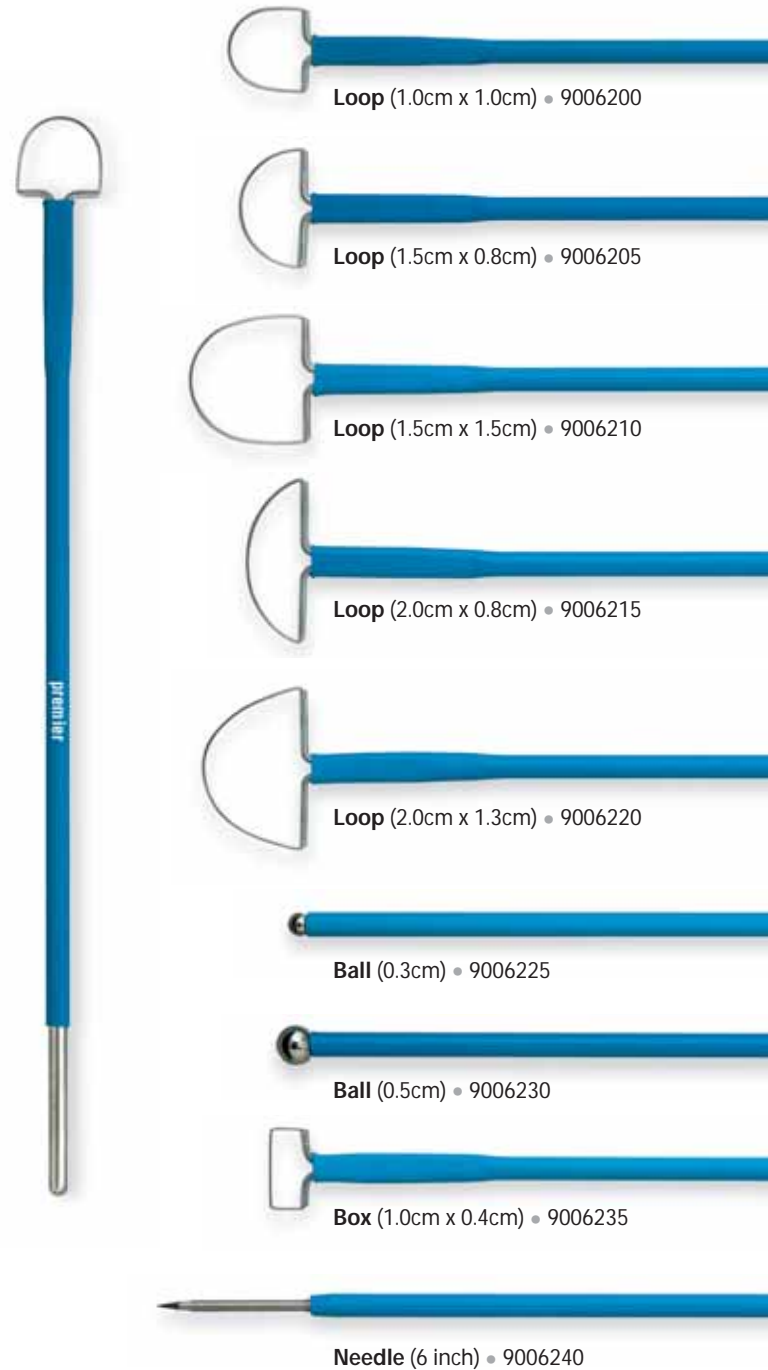

## "LOOP" Surgery Products

### "LOOP" Surgery Instruments

Designed specifically for use during LOOP and other electrosurgical procedures<sup>4</sup>, the non-conductive coating on these instruments ensures the safety of both patient and physician.

- Scratch Resistant, Non-conductive Coating
- Fully Autoclavable

### Use and Care Instructions

- Can be autoclaved or cold soaked
- Do not use with ultrasonic cleaner
- Do not use if coating is compromised
- Do not use with dry heat sterilization

### Coated Graves Specula

All Premier "Loop" Surgery specula are affixed with a smoke-removal tube. The tube, mounted on the upper blade, provides safe and convenient evacuation<sup>5</sup> of electrosurgical plume from the vaginal cavity.

9086165 Wide Medium (1 1/4" x 4 1/2" blades)

9086166 Wide Large (1 3/8" x 5" blades)

9086500 Small (3/4" x 3" blades)

9086501 Medium (1 1/4" x 4 1/2" blades)

9086502 Large (1 3/8" x 5" blades)



Wide Large  
(not actual size)

### Coated Lateral Vaginal Retractors



With Offset Blades  
(not actual size)

9086186 Lateral Retractor with offset blades

Fully opened, this retractor provides an additional 1" of retraction.

Overall Length - 10"  
Blade Size - 1" x 2 1/2"  
Fully opened dilation - 3 1/2"

9086182 Lateral Retractor (not shown)


Overall Length - 10"  
Blade Size - 1" x 2 1/2"  
Fully opened dilation - 2 1/2"

### Coated Hooks

(not actual size)

9086515 Iris Hook - 10"

9086520 Tenaculum Hook - 9 1/2"

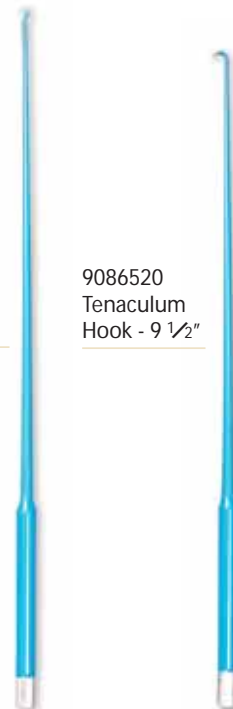

### Coated Tenaculum

(not actual size)



9086525 Schroeder Tenaculum - 9 3/4"

4. These instruments are not recommended for use with surgical lasers.  
5. Requires use of an air filtration system.

## Medicaments

### Monse's Solution



This thickened formulation of the popular hemostatic agent can be applied directly to the surgical site following excision of diseased tissue.



9045055 Box of twelve (12) 8ml bottles with applicators

### Lugol's Solution



Lugol's is a strong iodine contrast medium for identification of diseased cervical tissue.



9045056 Box of twelve (12) 8ml bottles with applicators

### Popular in LOOP, Biopsy and Colposcopic Procedures

Thickened formulation "adheres" to the surgical site for effective control of post-op bleeding.



## Pessary Products

### Pessary Products

For non-surgical management of uterine prolapse, cystocele and rectocele

- 100% silicone
- Easy to insert, remove and clean
- Single-patient use

**Ring Pessary** - This very common pessary most often is used for first- to second-degree prolapse. Ring Pessaries, available in nine sizes, are configured without and with support. Ring Pessaries with support also can be used on an accompanying cystocele. Application is accomplished easily simply by folding the ring before insertion.



Size	Diameter (inches)	Ring w/o Support Item #	Ring w/ Support Item #
#1	2.00	1040111	1040101
#2	2.25	1040112	1040102
#3	2.50	1040113	1040103
#4	2.75	1040114	1040104
#5	3.00	1040115	1040105
#6	3.25	1040116	1040106
#7	3.50	1040117	1040107
#8	3.75	1040118	1040108
#9	4.00	1040119	1040109

**Cube Pessary** - Designed for second- or third-degree prolapse or procidentia, as well as a cystocele and rectocele. Cube Pessaries feature drainage holes. Additionally, a silicone tie affixed to each cube aids in removal. Cube Pessaries are available in 10 sizes. Sizes range from 2.5cm (1") O.D. to 7.5cm (3") O.D.



Size	Width (cm)	Cube w/ Holes Item #
#0	2.5	1040400
#1	2.9	1040401
#2	3.2	1040402
#3	3.7	1040403
#4	4.1	1040404
#5	4.5	1040405
#6	5.5	1040406
#7	6.5	1040407
#8	7.0	1040408
#9	7.5	1040409

**Gellhorn Pessary** - Used for a second- to third-degree prolapse or procidentia. With drainage holes in the base, Gellhorn Pessaries are available in nine sizes. The knob of the gellhorn easily folds over for insertion and when properly inserted, rests on the posterior vaginal wall to provide support.



Size	Diameter (inches)	Gellhorn Item #
#0	1.50	1040200
#1	1.75	1040201
#2	2.00	1040202
#3	2.25	1040203
#4	2.50	1040204
#5	2.75	1040205
#6	3.00	1040206
#7	3.25	1040207
#8	3.50	1040208

**Donut Pessary** - For use with third-degree prolapse or procidentia as well as cystocele and rectocele. Donut Pessaries, while having a substantial amount of material between the inside and outside diameter, are made of soft material that is compressed easily to aid in insertion. The donut is available in seven different sizes.



Size	Diameter (inches)	Donut Item #
#0	2.00	1040300
#1	2.25	1040301
#2	2.50	1040302
#3	2.75	1040303
#4	3.00	1040304
#5	3.25	1040305
#6	3.50	1040306




**1040600 Pessary Fitting Set** - Set of six different rings (#2 through #7) used in pessary size selection for individual patients. The Fitting Set comes with a cross-reference chart that translates ring size into corresponding size of other pessary styles.

## Needle Adapters

### Reusable Stainless Steel Needle Adapters

These needle adapters allow for economical extension of an ordinary hypodermic syringe. The now "extended" needle easily reaches the paracervical and uterosacral areas to provide improved control for local infiltration of anesthetic agents. Slightly malleable for additional convenience. Four different lengths. Luer lock connectors.

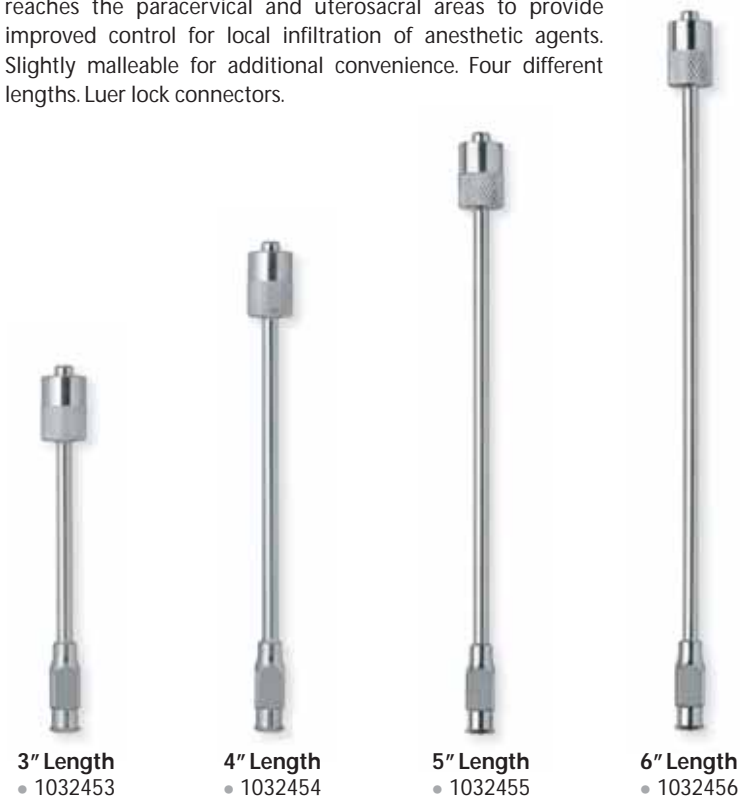


## German Biopsy Punch Forceps

German Precision Instruments

- Exact Bite for Clean Tissue Extraction
- Firm Hold of Biopsy Sample for Confident Procedures
- Spring Grip for Optimal Handling
- Swivel Lock to Ensure Sample Retention
- Popular 10" Length
- Outstanding Value, One-Year Unconditional Warranty

Made in Germany from high-grade surgical stainless steel, Premier®'s Biopsy Punch Forceps offer optimal performance and ultimate quality. Through superior craftsmanship, these four classic shapes afford the highest level of clinical care.

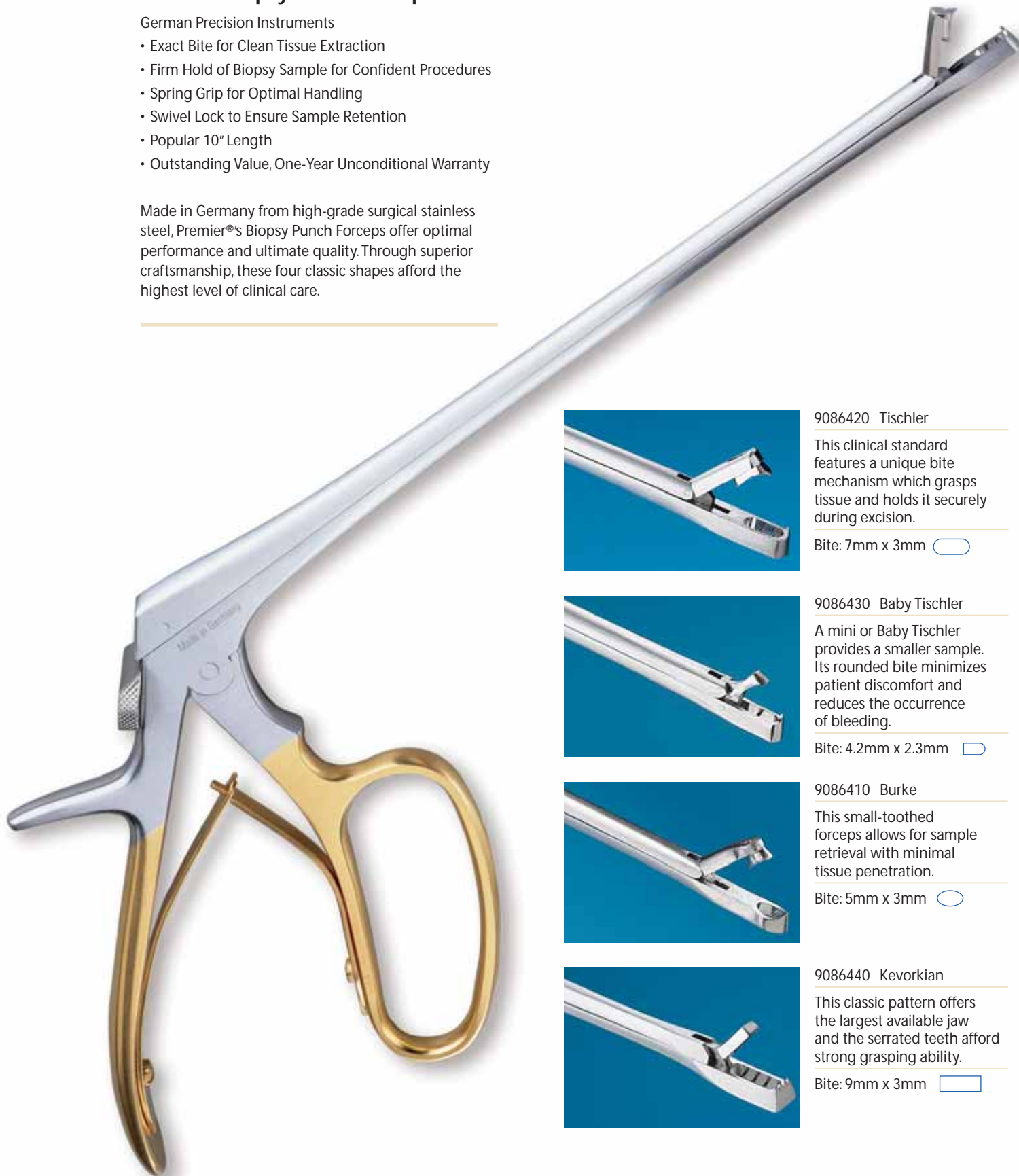

9086420 Tischler

This clinical standard features a unique bite mechanism which grasps tissue and holds it securely during excision.

Bite: 7mm x 3mm



9086430 Baby Tischler

A mini or Baby Tischler provides a smaller sample. Its rounded bite minimizes patient discomfort and reduces the occurrence of bleeding.

Bite: 4.2mm x 2.3mm



9086410 Burke

This small-toothed forceps allows for sample retrieval with minimal tissue penetration.

Bite: 5mm x 3mm



9086440 Kevorkian

This classic pattern offers the largest available jaw and the serrated teeth afford strong grasping ability.

Bite: 9mm x 3mm

# Specula



Graves Specula  
Upper Blade



Pederson Specula  
Upper Blade



## Graves

9086100 Small ( $\frac{3}{4}$ " x 3" blades)

9086101 Medium (1  $\frac{1}{4}$ " x 4  $\frac{1}{2}$ " blades)

Medium  
(not actual size)



## Pederson

9086110 Small ( $\frac{3}{4}$ " x 3" blades)

9086111 Medium (1  $\frac{1}{4}$ " x 4  $\frac{1}{2}$ " blades)

9086112 Large (1  $\frac{3}{8}$ " x 5" blades)

Medium  
(not actual size)

## Specula



### Weisman-Graves

9086140 LLO (1 1/4" x 4 1/2" blades)

9086141 RLO (1 1/4" x 4 1/2" blades)

LLO (not actual size)



Weisman-Graves  
Specula  
Upper Blade


### Diameters, French Scale

	mm	inch
○ 9	3.0	.118
○ 10	3.33	.131
○ 11	3.67	.144
○ 12	4.0	.158
○ 13	4.33	.170
○ 14	4.67	.184
○ 15	5.0	.197
○ 16	5.33	.210
○ 17	5.67	.223
○ 18	6.0	.236
○ 19	6.33	.249
○ 20	6.67	.263
○ 21	7.0	.275
○ 23	7.67	.302
○ 25	8.33	.328
○ 27	9.0	.354
○ 29	9.67	.380
○ 31	10.33	.407
○ 33	11.0	.433

## Uterine Dilators

**Hank** French scale (not actual size)

**Set of Six** • 1030541

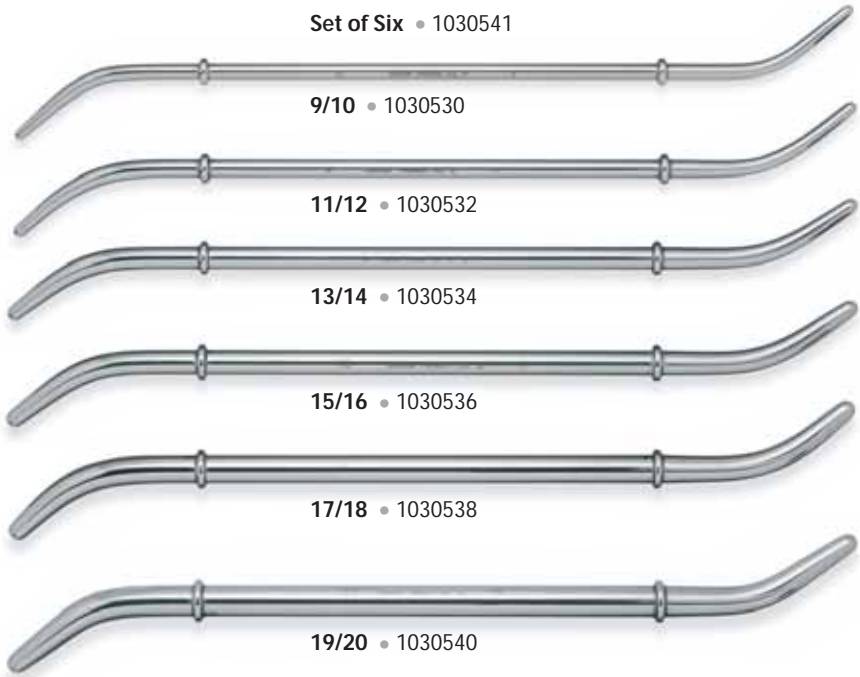

**9/10** • 1030530

**11/12** • 1030532

**13/14** • 1030534

**15/16** • 1030536

**17/18** • 1030538

**19/20** • 1030540

**Pratt** French scale (not actual size)

**Set of Eight** • 1030580



**13/15** • 1030565

**17/19** • 1030567

**21/23** • 1030569

**25/27** • 1030571

**29/31** • 1030573

**33/35** • 1030575

**37/39** • 1030577

**41/43** • 1030579

French	35	37	39	41	43
mm	11.67	12.33	13.0	13.67	14.33
inch	.459	.485	.512	.538	.564

## Uterine Dilators

**Hegar** (not actual size) **Set of Eight** • 9030559



3/4 (3mm/4mm) • 9030544



5/6 (5mm/6mm) • 9030546



7/8 (7mm/8mm) • 9030548



9/10 (9mm/10mm) • 9030550



11/12 (11mm/12mm) • 9030552



13/14 (13mm/14mm) • 9030554




15/16 (15mm/16mm) • 9030556



17/18 (17mm/18mm) • 9030558

**Kleegman** (not actual size)

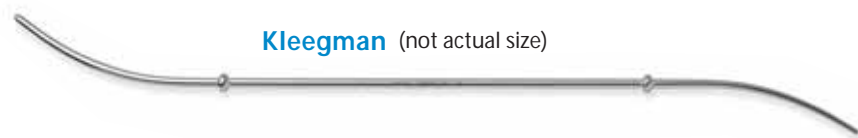

**Stainless** 1.8mm/2.5mm • 1030581



**Teflon** 1.5mm/2.1mm • 1030582

## Diameters, Millimeters

◦	1		11
○	2		12
○	3		13
○	4		14
○	5		15
○	6		16
○	7		17
○	8		18
○	9		
○	10		

## Cervical Dilators

Each Premier Os dilator is manufactured to exacting standards. Made of Teflon®, which is malleable, 'warm' and self-lubricating, these instruments deliver excellent patient comfort. Usage and care instructions are provided below.

**Os Locator** is distinguished with a single ring and a yellow band around the handle. It is particularly useful in locating the external Os in patients presenting cervical stenosis.


**Canal Finder** has a double ring and a blue band around the handle. The slightly longer tip insinuates the cervical canal, while producing dilation.

**Fundus Dilator** is fitted with a triple ring and a green band on its handle. Having the longest tip of the set, this particular instrument can be passed through the entire length of the cervical canal and will open the internal Os.

**The following guidelines ensure proper care and performance of the Os Dilators:**

- Do not bend the instrument more than 90° from its longitudinal axis.
- Inspect instrument for cracks or breakage prior to EACH use.
- Clean and sterilize before each use.
- Os Dilators can be cleaned and sterilized by any approved method. (Dry heat sterilization not recommended.)

**OS Dilators** Reusable (not actual size) **Set of Three** • 1030585

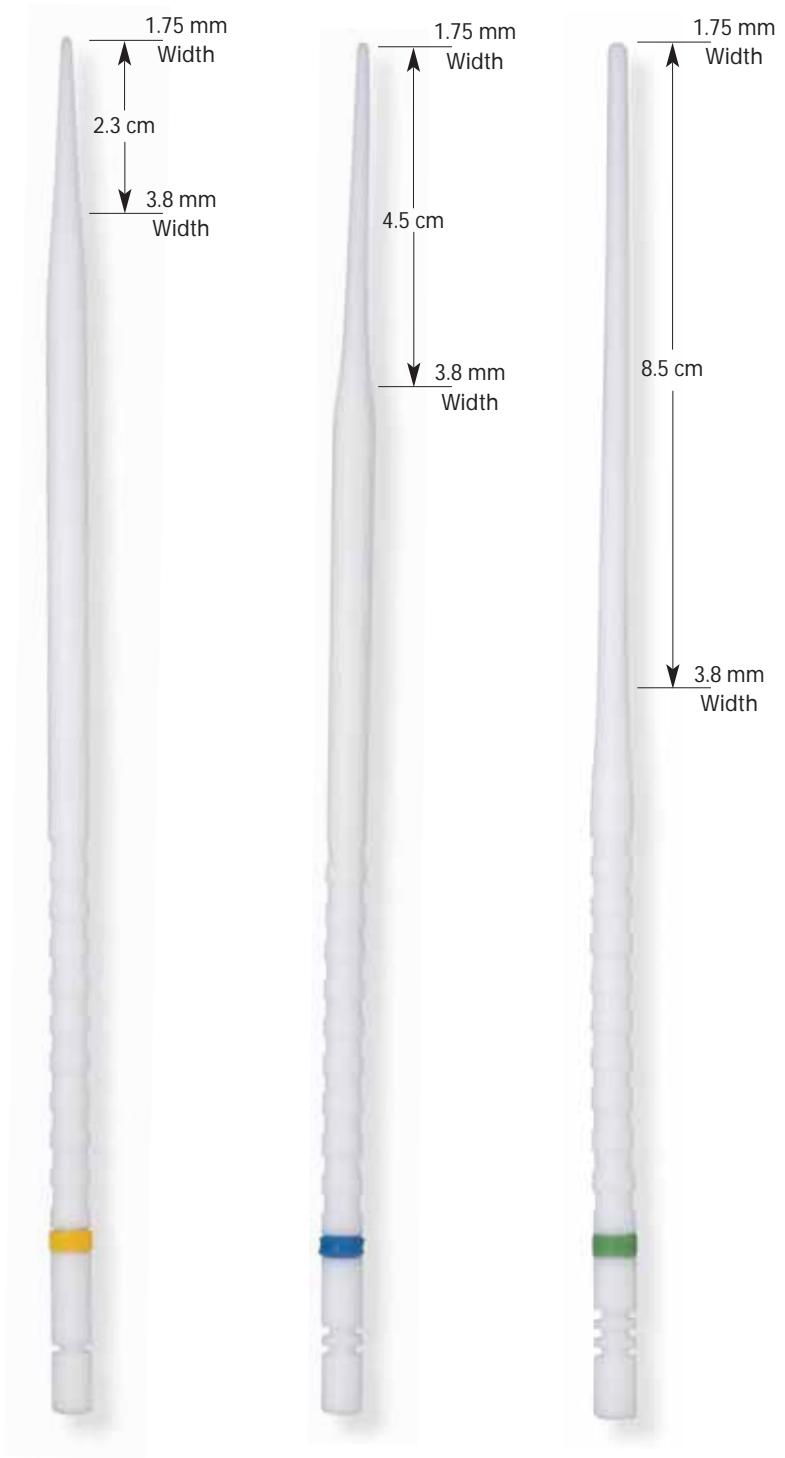

**OS Locator**  
• 1030586  
(Yellow Ring)

**Canal Finder**  
• 1030587  
(Blue Ring)

**Fundus Dilator**  
• 1030588  
(Green Ring)

## Uterine Sounds

### Sims



Graduated/cm • 1030522



Graduated/cm (2mm) • 1030524



Sim's Handle Style (not actual size)

All tips shown actual size.

All Premier Sounds are silver plated and malleable.

2mm (not actual size) All tips shown 200% actual size.

## Endometrial Biopsy Curettes

### Novak



2mm • 1030312

3mm • 1030313

4mm • 1030314



### Kevorkian-Young

(not actual size)



Malleable • 1030330  
(not shown)

Malleable with Basket • 1030331  
(not shown)



Rigid • 1030332



Rigid with Basket • 1030333

### Duncan


(not actual size)




Rigid • 1030343

# Instruments


## Atraumatic Jaws

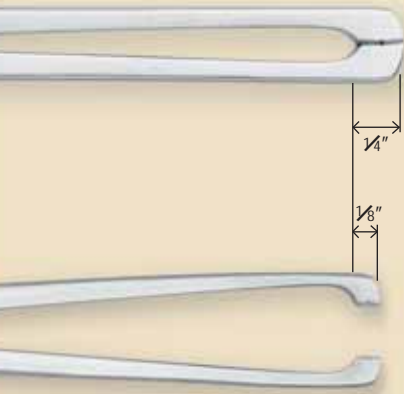

## Tenaculum (not actual size)

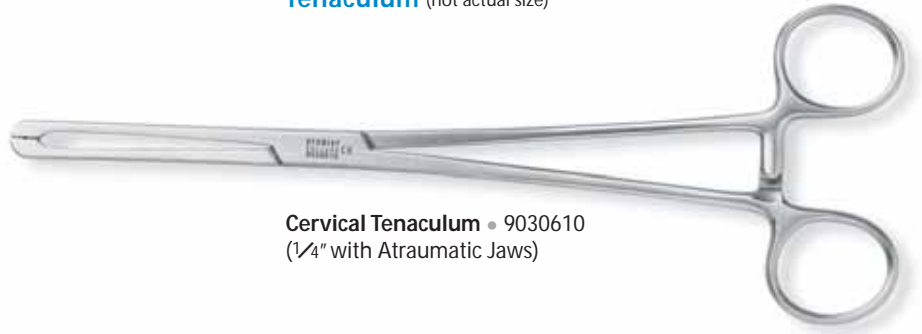


**Cervical Tenaculum • 9030610**  
(1/4" with Atraumatic Jaws)

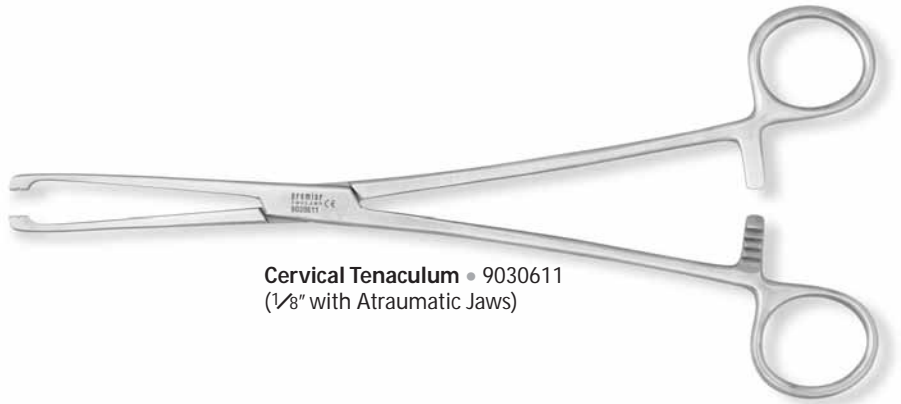


**Cervical Tenaculum • 9030611**  
(1/8" with Atraumatic Jaws)



**Emmett Tenaculum Hook #1**  
• 1030603 (9", Half Curve)

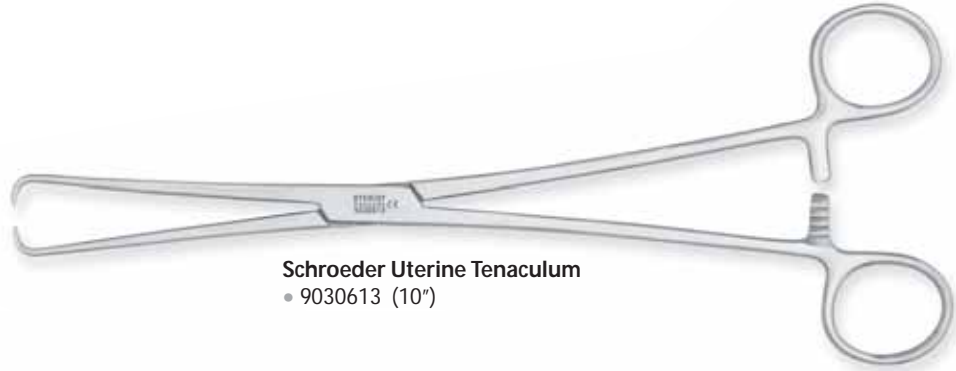

**Schroeder Uterine Tenaculum**  
• 9030613 (10")

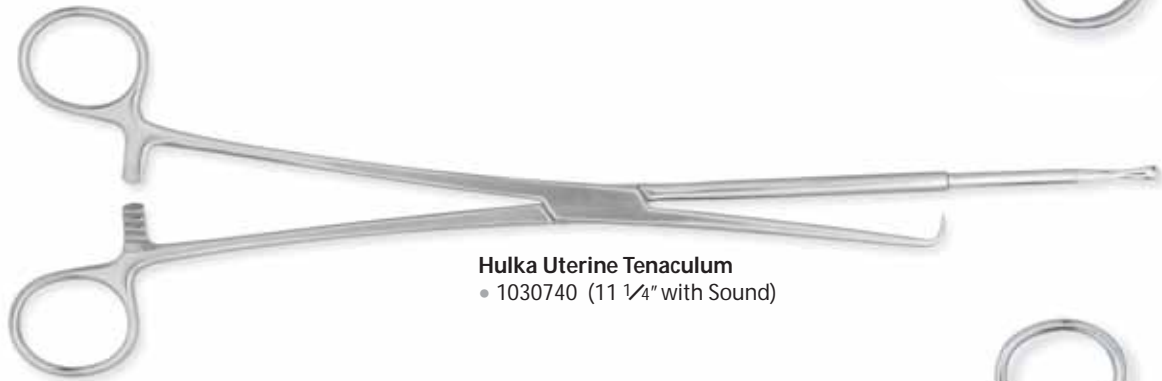

**Hulka Uterine Tenaculum**  
• 1030740 (11 1/4" with Sound)



**Hulka Uterine Tenaculum**  
• 1030742 (11 1/4" with Insulated Sound)



**Hulka Uterine Tenaculum**  
• 1030744 (11 1/4" with Atraumatic Jaw and Sound)

## Instruments

### Cannula Instruments

**Jarcho Self-Retaining Cannula** • 1030840  
(not actual size)



**Kahn Uterine Trigger Cannula**  
• 1030842 (with 3 Spare Tips)  
(not actual size)



These instruments are particularly useful in the administration of radiographic dyes and other solutions into the uterine cavity. Flexible malleable tips prevent the need for sounding or dilation.

### Replacement Parts for Cannula Instruments

1032100 Rubber Acorn (3/pack)

1032219 Teflon Tip, Short (3/pack)

1032220 Teflon Tip (3/pack)


1032221 Teflon Acorn

1032222 Collar/Screw (with Rubber Gasket)

1032223 Teflon Acorn for Corson

1032224 Collar Screw for Corson

Miscellaneous Instruments

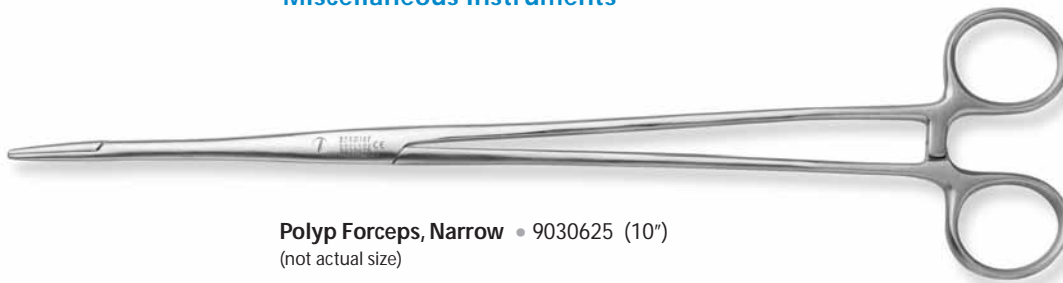

**Polyp Forceps, Narrow** • 9030625 (10")  
(not actual size)



**IUD Removal Hook** • 1030655  
(not actual size)



**IUD Removal Hook, Universal** • 1030656  
(not actual size)

Sheet Holders - OB/GYN

**Bethea Sheet Holder**  
• 1050222 (1 Pair)  
(not actual size)

